

WHEELCHAIRS





Everest & Jennings[®]

by graham-field

Value:

Graham-Field is known for a broad product offering at affordable pricing, without sacrifice to quality. Like you, we understand that the cost of a product is more than a price tag. The cost also includes the time and effort required to sell to and service your customers. At Graham-Field we are proud to offer the value found in the Everest & Jennings[®] line of wheelchairs.

Quality:

The E&J[™] name has been synonymous with quality manual folding wheelchairs since 1933. In a time when many manufacturers are producing products at the lowest possible cost, often cutting corners on quality, Graham-Field is standing by its dedication to quality through product design, specification control, testing and inspection. It is easy to see why the E&J Traveler[®] name is synonymous with quality.

Service:

Our products are often only a day away with centralized distribution located in the metro areas of Atlanta, Georgia; St. Louis, Missouri; and Los Angeles, California. Graham-Field also carries a full line of E&J parts and upholstery to meet your service needs. In addition, our customer service department and technical service teams are always ready to help with your product needs.



Transport Chairs

Aluminum Transport Chair

EJ78X-1 & EJ79X-1 Series

The Everest & Jennings Transport Chair is designed as an affordable short-term patient transfer. Its compact size and lighter weight allow a companion to lift and store the chair with ease.

- 17 and 19 inch (43.18 and 48.26 cm) seat widths [19 inch on Deluxe model only]
- Lightweight aluminum frame available in Aqua, Black, Red, Blue, Silver and Pink [Black and Blue Celeste are available for Deluxe Models only]
- Rear mounted attendant locks on rear casters make it easy for a companion to engage
- Fold-down back makes transport and storage easy
- Black, padded, nylon upholstery is flame retardant
- Tool-free adjustable detachable swingaway footrests
- Comfortable fixed full length arms with padded armrests
- Comes standard with 3 inch (7.62 cm) leg adjustable loops and strap to slide over installed footrests
- Includes carry straps on seat for lifting and folding ease
- 42 inch (106.68 cm) length positioning belt
- Weighs 20 lb (9.07 kg) [without front rigging]
- 250 lb (113.40 kg) weight capacity- EVENLY DISTRIBUTED
- One (1) year limited warranty

colorSelect



E1038*

Steel Transport Chair

EJ79X-1 Series

- Steel frame is finished in a durable hammertone coating
- Weighs 24 lb (10.89 kg) [without front rigging]
- 250 lb (113.40 kg) weight capacity - EVENLY DISTRIBUTED
- One (1) year limited warranty

E1038*



Deluxe Transport Chair

EJ87X-1 Series

- Handle mounted attendant locks
- 12 inch (55.88 kg) flat-free rear wheels with 8 inch (20.32 kg) swivel front casters
- Weighs 21.38 lb (9.70 kg) [without front rigging]
- 300 lb (136.08 kg) weight capacity - EVENLY DISTRIBUTED
- One (1) year limited warranty



Bariatric Transport

This chair satisfies additional compact transport needs with extra capacity. Handle-mounted attendant wheel locks are easily engaged by a companion.

- Dual crossbrace on gusseted aluminum frame with silver epoxy painted finish
- Available in 24 inch (60.96 cm) width, 18 inch (45.72 cm) seat depth
- 12 inch (30.48 cm) flat-free rear wheels with 8 inch (20.32 cm) swivel front casters
- Removable flip-back desk length padded arms provide user comfort
- Detachable swingaway footrests
- Available with anti-tippers on models after serial number IOS291
- Weighs 31 lb (14.06 kg) [without front rigging]
- 400 lb (181.44 kg) weight capacity - EVENLY DISTRIBUTED
- Three (3) year limited warranty
- Item No. EJ777-3



E1039*

Manual Folding Wheelchairs

Advantage® LX

3H02 Series

The Advantage, designed as a low-maintenance wheelchair, provides comfort and durability at an affordable price.

- Available in 16, 18 and 20 inch (40.64, 45.72 and 50.80 cm) seat widths, 16 inch (40.64 cm) seat depth
- Black padded nylon and vinyl upholstery with standard chart pocket on back
- Carbon steel frame with hammertone finish provides strength and durability
- Tool-free adjustable front rigging with standard heel loops on swingaway footrests
- Comfortable desk and full length arms with padded armrests
- 300 lb (136.08 kg) weight capacity - EVENLY DISTRIBUTED
- One (1) year limited warranty



K0001*

Traveler® SE Plus

3E05 Series

The Traveler SE Plus is the wheelchair of choice for standard transportation and hemi position.

- Black vinyl upholstery with standard chart pocket on back
- Two-position axle and caster fork provide 2 inch (5.08 cm) seat height adjustment
- 300 lb (136.08 kg) weight capacity - EVENLY DISTRIBUTED
- One (1) year limited warranty



K0001*
K0002

Traveler® LX

3E02 Series

The Traveler LX incorporates the legendary Everest & Jennings wheelchair chrome finish frame with saddle-joint tubular construction with dual axle. The premium blue vinyl upholstery is reinforced with an extra foam and heavy duty inner liner for enhanced comfort and durability.

- Chrome finish provides a long lasting and easy to clean finish
- Gray non-marking tires and aluminum footplates
- 300 lb (136.08 kg) weight capacity - EVENLY DISTRIBUTED
- Two (2) year limited warranty



K0001*
K0002

Lightweight Folding Wheelchairs

Traveler® L3 Plus and Traveler® L3 Plus XVI

3F012, 3F013 & 3F014 Series

The Traveler L3 Plus is an improved version of the previous Traveler L3 series. The wheelchair comes standard with adjustable anti-tippers and a fully detachable 2 inch (5.08 cm) seat extension with hardware.

- Available in 16, 18 and 20 inch (40.64, 45.72 and 50.80 cm) seat widths including 2 inch (5.08 cm) removable seat extension kit converting the chair from 16 inch depth to 18 inch depth if necessary
- Two position axle and caster fork provides 18 to 20 inch seat height in regular models and seat height 16.30 to 17.90 inch (41.40 to 45.47 cm) in the XVI models
- Black nylon upholstery with standard chart pocket on back
- Includes fastener tab under seat connecting both seat and back upholstery to eliminate gap
- Non-marking polyurethane tires with 24 inch (60.96 cm) rear wheels with 8 inch (20.32 cm) casters in the regular models and 20 inch rear wheels with 6 inch (15.24 cm) casters in the XVI models
- Anti-tippers come standard
- Flip-back padded arms available in fixed desk length, adjustable height desk and detachable flip-back full length
- Hammertone finish and carbon steel frame provide strength and durability
- Tool-free adjustable front rigging with heel loops standard on swingaway footrests
- 300 lb (136.08 kg) weight capacity - EVENLY DISTRIBUTED
- Five (5) year limited warranty



24 inch (60.96 cm) wheel model shown



24 inch (60.96 cm) wheel model shown



20 inch (50.80 cm) wheel model shown

Traveler® L4

3F02 & 3F03 Series

The Traveler L4 is a high-strength lightweight chair that weighs less than 34 lb (15.42 kg) [without front rigging] with an adjustable back.

- Available in 16, 18 and 20 inch (40.64, 45.72 and 50.80 cm) seat widths and 16 inch seat depth. Plus 18 and 20 inch seat widths, 18 inch seat depth
- Quick release axle and height adjustable desk arms are available for some models
- Black nylon upholstery with standard chart pocket on back
- Hammertone paint provides a long lasting and durable finish
- Anti-tippers come standard
- Comfortable flip-back desk and full length arms with padded armrests
- 250 lb (113.40 kg) weight capacity - EVENLY DISTRIBUTED
- 300 lb (136.08 kg) weight capacity - EVENLY DISTRIBUTED on models with 18 inch seat depth
- Five (5) year limited warranty



K0004*

Bariatric Folding Wheelchairs

Traveler® HD

3G01 Series

The Traveler HD is a heavy duty wheelchair featuring a reinforced frame, double crossbrace, and a long lasting, embossed back and seat vinyl upholstery with a double inner liner for durability and support.

- Available in 20, 22 and 24 inch (50.80, 55.88 and 60.96 cm) seat widths, 18 inch (45.72 cm) seat depth (1 inch / 2.54 cm longer than standard)
- Embossed black vinyl upholstery with standard chart pocket on back
- Hammertone paint provides a long lasting and durable finish
- Reinforced frame at all critical load-bearing points
- Two-position axle and caster fork provides 2 inch (5.08 cm) seat height adjustment
- Tool-free adjustable front rigging with standard heel loops on footrests and durable aluminum black footplates
- Comfortable full and desk length arms with padded armrests
- 500 lb weight capacity - EVENLY DISTRIBUTED
- Five (5) year limited warranty



Traveler® HTC

3J01 Series

The Traveler HTC (Hospital Transport Chair) is an all-in-one solution wheelchair for the clinical transport need, as it incorporates a special attendant push handle and rear mounting attendant wheel locks for better caregiver assistance.

- Fixed desk arms and front rigging lock nuts ensure inventory control
- Anti-tippers and one O2/IV telescopic holder come standard
- 500 lb (226.80 kg) weight capacity - EVENLY DISTRIBUTED
- Five (5) year limited warranty



Paramount XD

5PX10 Series

The Paramount XD, built for the bariatric wheelchair user, is constructed of 1 inch (2.54 cm) tubing to meet high weight capacity demands.

- Available in 26 and 30 inch (66.04 and 76.20 cm) seat widths, 19 ½ inch (49.53 cm) seat depth
- Comfortable desk arms with padded armrests
- Front rigging with extra wide aluminum footplate
- 8 x 2 inch (20.32 x 5.08 cm) non-marking solid casters and 24 inch (60.96 cm) solid rear wheels with chromed handrims
- Black reinforced nylon upholstery with standard chart pocket on back
- 650 lb (294.84 kg) weight capacity - EVENLY DISTRIBUTED
- Five (5) year limited warranty



Specialty Folding Wheelchairs

Advantage® Recliner

3K01 & 3K02 Series

The Advantage Recliner, with a higher back and removable, full-support headrest for additional neck support, adjusts easily to fit every need.

- Available in 16, 18 and 20 inch (40.64, 45.72 and 50.80 cm) seat widths, 17 inch (43.18 cm) seat depth and in 22 inch (55.88 cm) seat width, 19 inch (55.88 cm) seat depth
- Durable hammertone frame features manual full reclining back range from 90-170 degrees
- Comfortable desk and full length arms with padded armrests
- Metal spreader bar on back which works as anti-fold device
- Embossed black vinyl upholstery with standard chart pocket on back
- Anti-tippers standard
- Tool-free adjustable elevating front legrests
- 450 lb (204.12 kg) weight capacity on 22 inch wide models - EVENLY DISTRIBUTED
- 300 lb (136.08 kg) weight capacity on 16, 18 & 20 inch wide models - EVENLY DISTRIBUTED
- Five (5) year limited warranty



K0001*
E1226
E0966
E0971
K0195

K0007*
E1226
E0966
E0971
K0195
E2201

Rehab Shower Commode

12022 & 12023 Series

The versatile E&J Folding Rehab Shower Commode chair has earned the reputation of being the premier self-propelled chair for personal hygiene care. Designed with a corrosion-resistant stainless steel frame, the chair can be positioned over a standard toilet or used as a portable, self-contained commode.

- Available with durable 24 inch (60.96 cm) rear wheels with handrims and solid rubber tires
- Removable full length arms with composite armrests
- Available in Low-back 14 ¼ inch (36.83 cm) on 16 inch (40.64 cm) and Hi-back 17 ¼ inch (43.82 cm) on 19 ½ inch (48.26 cm)
- Available in 17 inch (43.18 cm) seat with 18 inch (45.72 cm) seat depth on Low-back 14 ¼ and 17 ¼ inch and the wider version 20 inch (50.80 cm) seat width on 16 inch Low-back and Hi-back 19 ½ inch models
- Sealed vinyl cushioned oval seat and frame side contours allow easy access and comfort in toileting activities and bowel management procedures
- Removable, hygienic plastic pail and lid are included
- Non-corrosive composite footplates
- 250 lb (113.40 kg) weight capacity -EVENLY DISTRIBUTED on 14 ¼ and 17 ¼ inch back models
- 300 lb (136.08 kg) weight capacity -EVENLY DISTRIBUTED on 16 and 19 ½ inch back models
- Two (2) year limited warranty



E0165*



Configurations & Accessories

Traveler®

		Seat Width		Seat Depth		Flip Back Desk	Flip Back Full	Fixed Full	Adjust. Ht. Desk	Swing Foot	Elevating Leg
		inch	cm	inch	cm						
3F012120	Traveler L3 Plus	18.00	45.72	16.00	40.64	■				■	
3F012130	Traveler L3 Plus	18.00	45.72	16.00	40.64	■					
3F012140	Traveler L3 Plus	18.00	45.72	16.00	40.64		■			■	
3F012220	Traveler L3 Plus	16.00	40.64	16.00	40.64	■					■
3F012320	Traveler L3 Plus	20.00	50.80	16.00	40.64	■					■
3F012340	Traveler L3 Plus	20.00	50.80	16.00	40.64		■			■	
3F013120	Traveler L3 Plus	18.00	45.72	16.00	40.64			■		■	
3F013320	Traveler L3 Plus	20.00	50.80	16.00	40.64			■		■	
3F014120	Traveler L3 Plus XVI	18.00	45.72	16.00	40.64	■				■	
3F014220	Traveler L3 Plus XVI	16.00	40.64	16.00	40.64	■				■	
3F014320	Traveler L3 Plus XVI	20.00	50.80	16.00	40.64	■				■	
3F020120	Traveler L4	18.00	45.72	16.00	40.64	■				■	
3F020150	Traveler L4	18.00	45.72	16.00	40.64		■				■
3F020140	Traveler L4	18.00	45.72	18.00	45.72		■				■
3G010320	Traveler HD	20.00	50.80	18.00	45.72	■					■
3G010340	Traveler HD	20.00	50.80	18.00	45.72		■			■	
3G010350	Traveler HD	20.00	50.80	18.00	45.72		■				■
3G010420	Traveler HD	22.00	55.88	18.00	45.72	■					■
3G010450	Traveler HD	22.00	55.88	18.00	45.72		■				■
3G010520	Traveler HD	24.00	60.96	18.00	45.72	■					■
3G010540	Traveler HD	24.00	60.96	18.00	45.72		■				■
3J010330	Traveler HTC	20.00	50.80	18.00	45.72	■					■
3J010340	Traveler HTC	22.00	55.88	18.00	45.72	■					■
3J010530	Traveler HTC	24.00	60.96	18.00	45.72	■					■
3E050120	Traveler SE Plus	18.00	45.72	16.00	40.64	■					■
3E050130	Traveler SE Plus	18.00	45.72	16.00	40.64	■					■
3E050220	Traveler SE Plus	16.00	40.64	16.00	40.64	■					■
3E050330	Traveler SE Plus	20.00	50.80	16.00	40.64	■					■

Paramount

		Seat Width		Seat Depth		Flip Back Desk	Flip Back Full	Fixed Full	Adjust. Ht. Desk	Swing Foot	Elevating Leg
		inch	cm	inch	cm						
5PX10620	Paramount XD	26.00	45.72	19.50	49.53	■					■
5PX10830	Paramount XD	30.00	50.80	19.50	49.53	■					■



Cup Holder

EJCUPH19 Cup Holder




Oxygen & IV

B10320F	O2/IV Combination Holder w/2 hooks Fits: Advantage, SE, L3, L3 Plus, L4, HD, LX, Paramount XD
B10325	O2 Holder Fits: Advantage, SE, L3, L3 Plus, L4, HD, LX, Paramount XD
B10321F	IV Pole Holder w/ 2 hooks Fits: Advantage, SE, SE Plus, L3, L3 Plus, HD, LX, Paramount XD
B103206	IV Pole Holder w/ 1 hook Fits: Advantage, LX, SE, SE Plus, L3, L3 Plus, HD, LX, Paramount XD
EJ02BAG07	O2 Bag



Footrest

90763030 / 90763030A	with Composite Footplate
90763430	with Aluminum Footplate

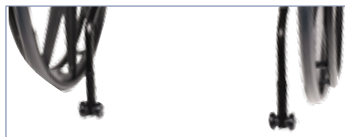

Anti-tippers

90763565	Fits: Advantage LX
90763065	Fits: Traveler SE & HD
90763065A	Fits: Traveler L3 Plus
90763233A	Fits: Traveler L3 Plus XVI

Advantage

		Seat Width		Seat Depth		Flip Back Desk	Flip Back Full	Fixed Full	Adjust. Ht. Desk	Swing Foot	Elevating Leg
		inch	cm	inch	cm						
3H020100	Advantage LX (NYLON)	18.00	45.72	16.00	40.64					■	
3H020120	Advantage LX (NYLON)	18.00	45.72	16.00	40.64	■					■
3H020230	Advantage LX (NYLON)	16.00	40.64	16.00	40.64	■					■
3H020140	Advantage LX (NYLON)	18.00	45.72	16.00	40.64		■			■	
3H021100	Advantage LX (VINYL)	18.00	45.72	16.00	40.64			■		■	
3H021120	Advantage LX (VINYL)	18.00	45.72	16.00	40.64	■					■
3H021130	Advantage LX (VINYL)	18.00	45.72	16.00	40.64	■					■
3H021320	Advantage LX (VINYL)	20.00	50.80	16.00	40.64	■					■
3K010130	Advantage Reclining	18.00	45.72	17.00	43.18	■					■
3K010150	Advantage Reclining	18.00	45.72	17.00	43.18		■				■
3K010230	Advantage Reclining	16.00	40.64	17.00	43.18	■					■
3K010350	Advantage Reclining	20.00	50.80	17.00	43.18		■				■
3K020430	Advantage Reclining	22.00	55.88	19.00	48.26	■					■
3K020450	Advantage Reclining	22.00	55.88	19.00	48.26		■				■

Transport Chairs

		Back Height	Seat Width	Seat Depth	Flip Back Desk	Flip Back Full	Fixed Full	Adjust. Ht. Desk	Swing Foot	Elevating Leg
EJ777-3	Bariatric Aluminum - Silver	24.00	60.96	18.00	45.72	■				■
EJ782-1	Aluminum - Red	19.00	48.26	16.00	40.64			■		■
EJ784-1	Aluminum - Silver	17.00	43.18	16.00	40.64			■		■
EJ786-1	Aluminum - Black	19.00	48.26	16.00	40.64			■		■
EJ788-1	Aluminum - Blue	19.00	48.26	16.00	40.64			■		■
EJ870-1	Deluxe Aluminum - Black	19.00	48.26	16.00	40.64			■		■
EJ796-1	Steel - Silver Vein	19.00	48.26	15.75	40.01			■		■


Rehab Shower Commode

		Back Height		Seat Width		Seat Depth		Arms	Rests
		inch	cm	inch	cm	inch	cm		
12022010	LoBack	14.25	36.20	17.00	43.18	18.00	45.72	■	■
12022110	HiBack	17.25	43.82	19.00	48.26	18.00	45.72	■	■
12023010	LoBack	16.00	40.64	19.00	48.26	16.00	40.64	■	■
12023110	HiBack	19.00	48.26	17.00	43.18	16.00	40.64	■	■

WARNINGS: Cancer and Reproductive Harm - www.P65Warnings.ca.gov

*For a complete list of sizes and options available, please visit www.grahamfield.com

All sales subject to Graham-Field terms and conditions (available at www.grahamfield.com)

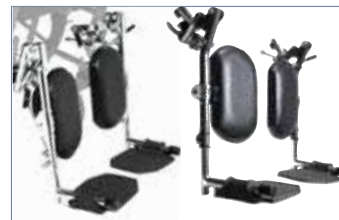

Legrest

90763031 / 90763031A	Elevating with Composite Footplate
90763431	Elevating with Aluminum Footplate
90763432	Articulating, Heavy-duty, with aluminum footplate & panel



Wheelchair Positioning Tray

EJTRAY18 Fits: All Models



Padded Seat

90835070	Oval cutout seat for 17 1/4 & 14 1/4 inch (44.45 & 36.83 cm)
90835071	Oval cutout seat for 16 & 19 inch (40.64 & 48.26 cm)
90835175	Open front (C-shaped) for 17 1/4 & 14 1/4 inch (44.45 & 36.83 cm)



Positioning Belts

300-505	with Auto Style
EJ76-BELT	with Plastic Buckle

Color variations may occur due to monitor or printer settings.

*1 770.368.4700

Information contained herein is subject to change without notice. The most current and complete product information can be found on our website.

www.grahamfield.com



© 2014, GF Health Products, Inc. All Rights Reserved. Everest & Jennings is a trademark of GF Health Products, Inc. GF Health Products, Inc. is an ISO 13485 Certified Company.

grahamfield

GF1400066Rev120